

Install Guide

Wall side first, then the bed.

BlackJack
by HATCH+MED


TOOLS Phillips screwdriver • Flat-head screwdriver

Scan for more info



Wall side


Choose one of these two ways to connect BlackJack to the 37-pin wall port.

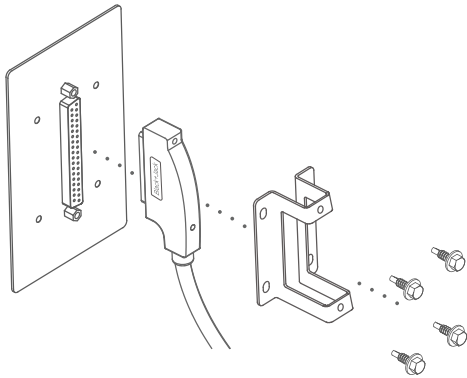

Right-angle wall cable CFR2 · CFR6

1. Plug the cable into the 37-pin port so the cable exits downward.
2. Tighten both thumbscrews.
3. Secure the cable with the metal bracket, a P-clamp about 1 foot below the port, or both.


or

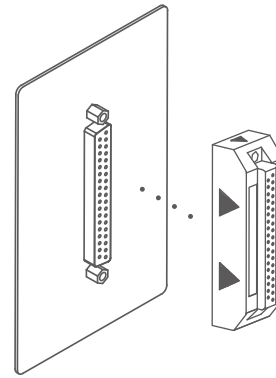

Magnetic wall adapter XA

1. Plug the adapter into the 37-pin port with the arrow on top.
2. Tighten both thumbscrews.

Bed side

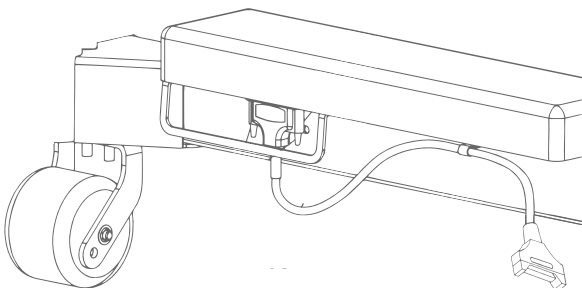

1. Plug the bed cable into the bed's 37-pin port, under or behind the bed.
2. Tighten both thumbscrews.
3. Anchor the cable to the bed frame with a P-clamp as a second hold.
4. Bring the magnetic ends together to connect.

Cleaning Instructions

Routine cleaning, plus a yearly check of each magnetic connector.

BlackJack
by HATCH+MED



BlackJack maintenance kit HMNO-4752

3 brushes, 2 tins of magnetic putty, 3 WD-40 precision pens and a microfiber cloth.



Scan for cleaning
and maintenance help

Clean the cable

1. Clean your hands and put on gloves.
2. Pre-clean the cable with a cloth moistened, not soaked, with mild detergent and water.
3. Wipe the cable with disinfectant on a cloth or microfiber cloth, or use a cleaning wipe.
4. Let the disinfectant dry, then wipe the cable with a cloth dampened with clean water.

APPROVED CLEANERS

Use one of these at the proper dilution to clean, deodorize or disinfect BlackJack cables.

Triple S Sanotracin RTU

JohnsonDiversey Alpha HP / MC

Virex II 256

Clorox Healthcare Bleach Germicidal Cleaner

Clorox Germicidal Wipes

Ecolab Asepti-wipe II

Triple S Perisept

IMPORTANT

Never put cleaning solution on the 37-pin connectors or the magnetic connectors. Fluid can get into the connection and damage the cable.

Yearly maintenance

Inspect and clean each magnetic connector, wall side and bed side, at least once a year to keep the factory warranty.

1. Disconnect the bed cable and clear the cord-out alarm. Replace any cable with damage, such as a cut jacket.
2. Lift debris off each magnetic connector with the putty, then clean every pin with the WD-40 pen until each shows clean gold.
3. Brush the connector with the nylon end of the brush. On the bed cable, also drag the rigid edge across the pins, five sweeps each way.
4. Wipe off extra WD-40 with the microfiber cloth and check the mounting screws.
5. Reconnect, test the nurse call and bed exit calls, and reset both.
6. Record the date and technician for your warranty records.