


iPad Assembly

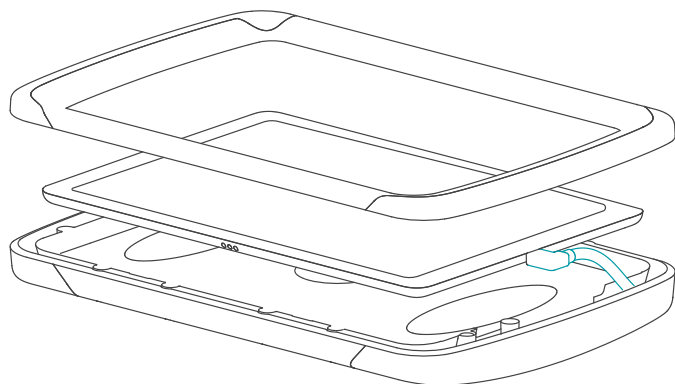

1. Connect the iPad to the charging cable and lay it flat.
2. Seal the front cover, flip the unit over, and secure the case screws, **do not insert the kickstand screws yet.**
3. After assembly, verify that the volume and power buttons work as expected.


Security Screws

HMNF-7507

Use a **1.5 mm / #4 spanner tool** and the right screw length for each hole for a flush, sealed finish. Using the **wrong screw length** can **damage** the unit.




Use provided
security screw
bit

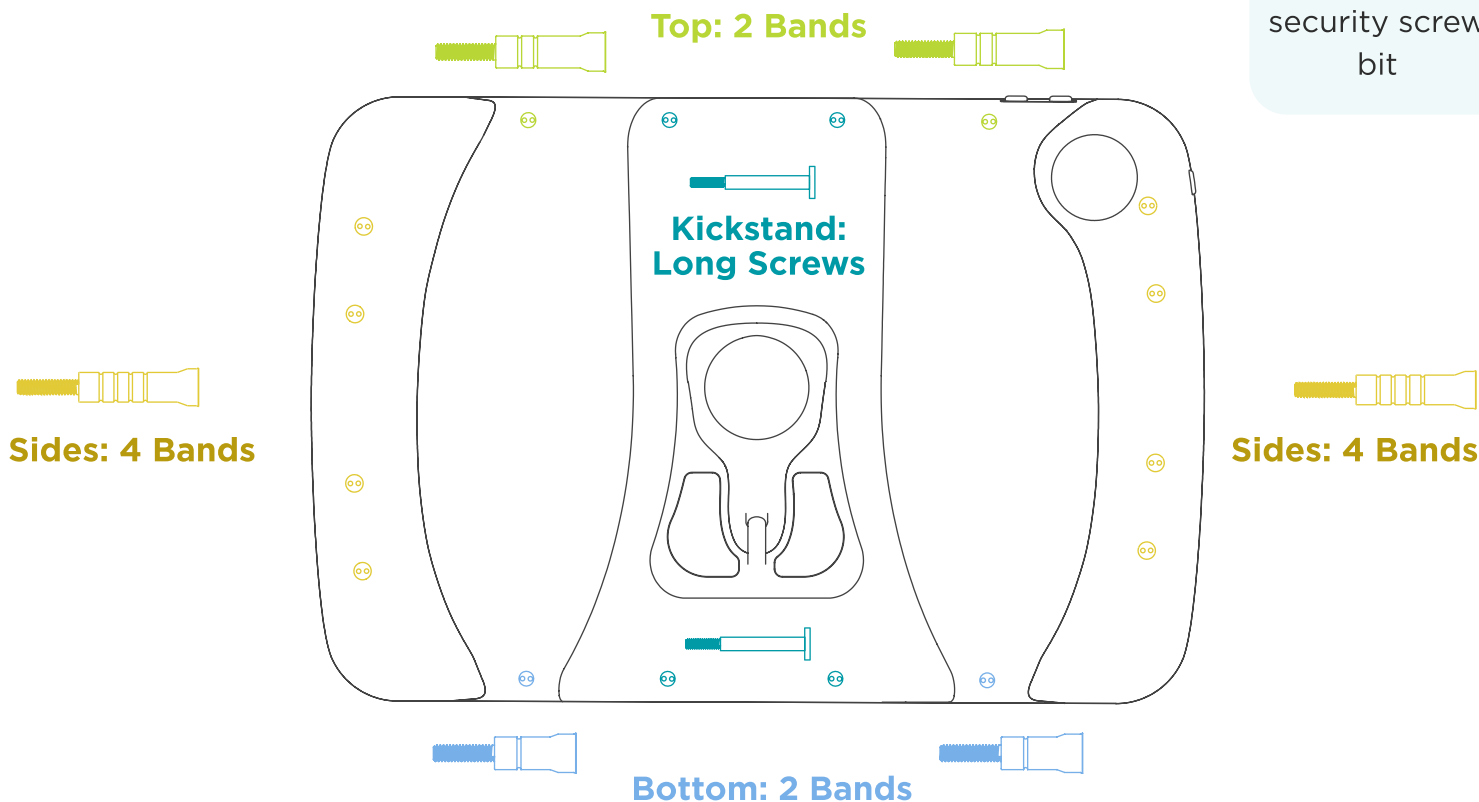

If using an impact driver, set your driver to its **lowest torque setting.**

Kickstand & Rail Clip

HMNF-7505

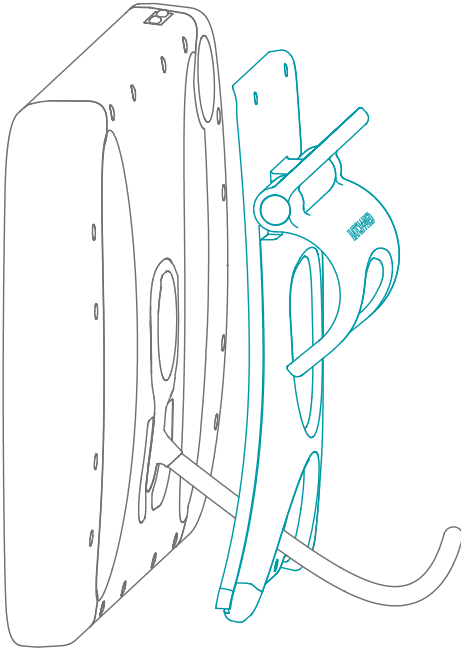


1. Align the **kickstand** screw holes with the matching holes on the back of the case.
2. Drive each kickstand screw into its matching hole using the band chart on the previous page.

Make sure to use the longer kickstand screws and set the torque to the lowest setting.

Nurse Call (Optional)


HMNF-7002A

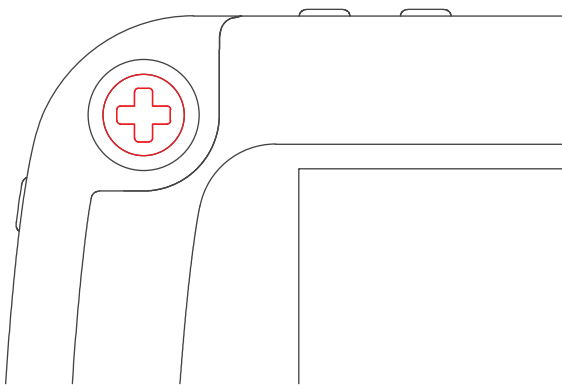

Certain PillowCase models include a physical UL 1069 **nurse call button** in the top-left corner. Pressing the red button or disconnecting the magnet places an **emergency call** when connected to HighLightX Smart Hub.

USB-C Port (Optional)

HMNF-7004B & HMNF-7004C

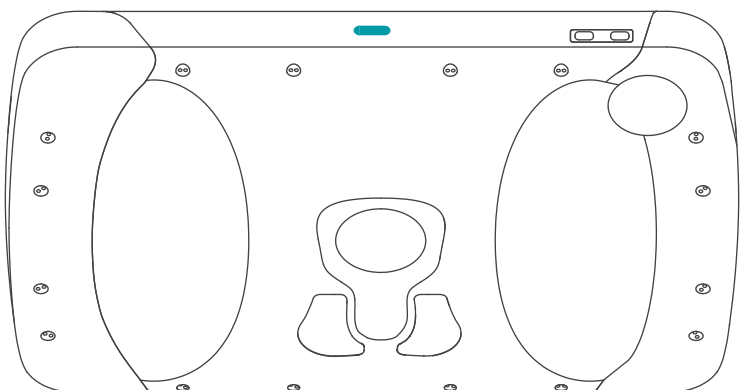


Certain PillowCase models include a waterproof USB-C port on the top edge for iPad MDM access or alternate charging - make sure to use only Apple-certified USB-C power supplies.

Power & Security Lock

Security Lock

HMNF-7509

Install the security lock on the **PillowCase** side of the cable, so the magnet can break away and slide up to the lock.

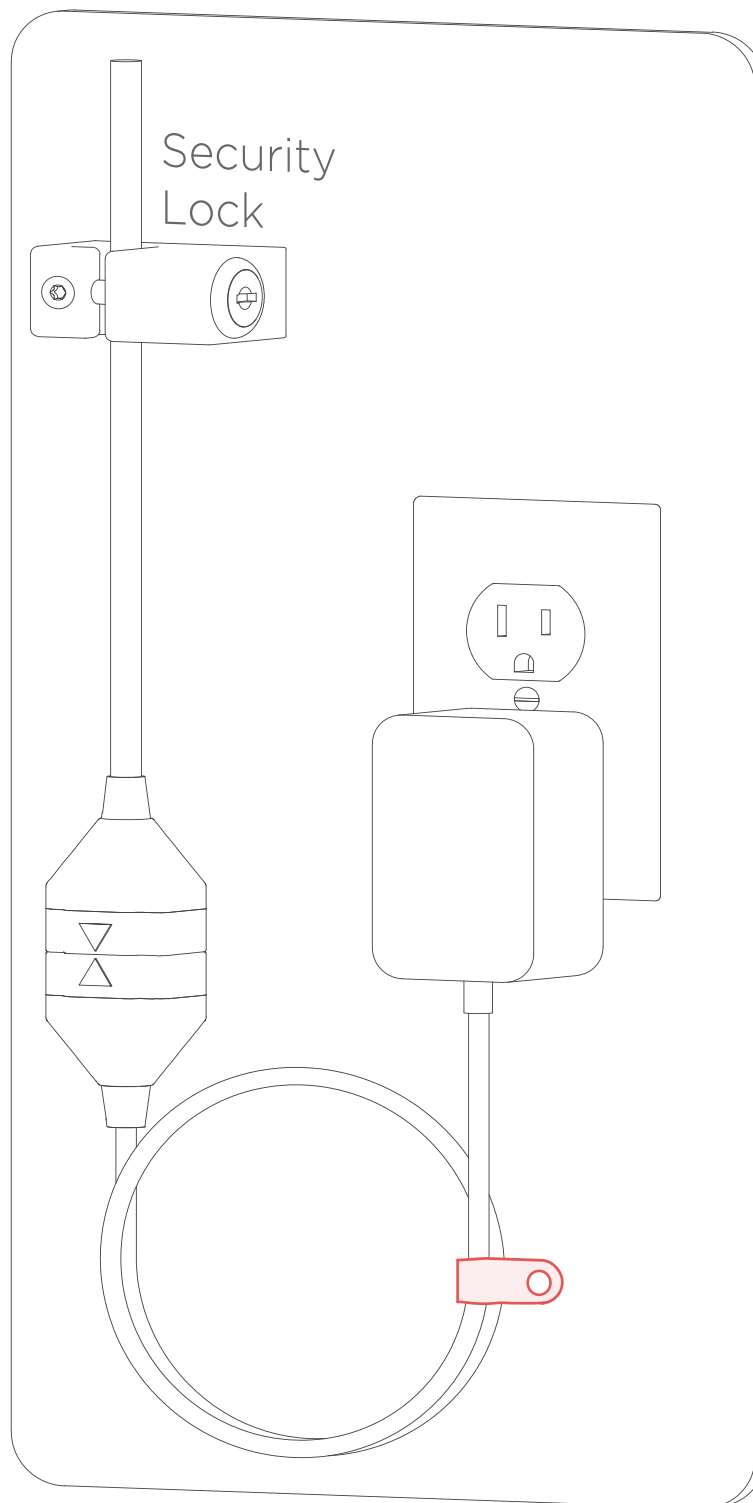

Power

Connect to the outlet. Confirm the iPad is charging and the outlet is not loose.

A **P-clamp** is required on the **power supply side** of the magnetic cable.

Cleaning Guide

Cleaning Steps

1. Clean your hands and wear gloves to avoid contact with solutions.
2. Use a liquid disinfectant or a cleaning wipe.
3. Pre-clean the cable with a cloth dampened (not soaked) in mild detergent and water.
4. Apply the proper dilution to a cloth or microfiber cloth, or use a wipe.
5. Disinfect the cable and case.
6. Allow to dry.
7. Wipe with a cloth soaked in clean water.

Approved Solutions

- Triple S Sanotracin RTU™
- JohnsonDiversey Alpha HP™/MC
- Virex II™ 256
- Clorox Healthcare® Bleach Germicidal Cleaner
- Clorox Germicidal Wipes
- Ecolab Asepti-wipe I
- Clorox™ or Diversey™ H₂O₂ Wipes

Do **not** clean the 14-pin magnet connectors with solution. To clear debris, pull the magnet apart and brush it with a nylon brush.

Contact our Success Team for installation questions, spare parts and warranty claims

844-428-2463 techsupport@hatchmed.com

