

PillowCase

by HATCH+MED



Transform your iPad into an FDA medical device - PillowCase keeps the iPad always powered, secured, and safe from fluid or chemical damage around the patient.




Compatible with iPad 7-10th Gen,
iPad Air 3-5th Gen
& iPad mini



Optional kickstand with a
spring loaded wall & bed rail
clip



UL listed plug & play
integration with any nurse call
system

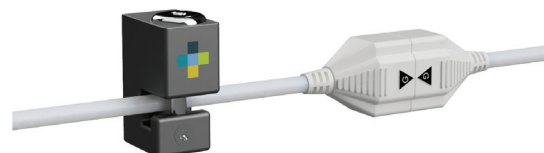

Wall mounted lock allows
magnet to disconnect while
preventing iPad theft



info@hatchmed.com



844-428-2463

PillowCase Configurations



Version A: iPad Security Case

- No physical nurse call buttons on case
- No software integration required
- Direct connection to power outlet with UL listed power supply
- Magnetic cable with 8ft on case side & 2ft on wall side
- Water-resistant and cleaning agent proof enclosure
- Kickstand & security lock accessories



Version B: Emergency Call Button

- UL 1069 hardware button for emergency calls on any nurse call system via 1/4" port
- No software integration required
- Wall mounted hub with UL listed nurse call connection and power supply
- Magnetic cable with 8ft on case side & 2ft on wall side
- Water-resistant and cleaning agent proof enclosure
- Kickstand & security lock accessories



Version C: Digital Pillow Speaker

- UL 1069 hardware button for emergency calls
- Software room controls, TV controls, clinical requests directly through nurse call system patient station (no API setup required) for Responder 5, HillRom and other systems
- Wall mounted hub with UL listed nurse call connection and power supply
- Magnetic cable with 8ft on case side & 2ft on wall side
- Water-resistant and cleaning agent proof enclosure
- Kickstand & security lock accessories



PillowCase Specifications

Size & Weight

| | Larger Case | Smaller Case |
|------------------------|--|---|
| Supported iPads | iPad 7th, 8th, 9th & 10th Gen iPad Air 3rd, 4th & 5th Gen iPad Pro 10.5" & 11" | iPad mini 6th Gen |
| Weight | 1.0 lb (without iPad) 2.1 lb (with iPad) | 0.8 lb (without iPad) 1.5 lb (with iPad) |
| Dimensions | 11.5" W x 7.9" H x 1.3" D | 9.7" W x 6.3" H x 1.3" D |

General Info

| | |
|----------------------------------|---|
| Connection & Charging | 10ft keyed magnetic cable (8ft on PillowCase side & 2ft on wall side) |
| Materials | Polycarbonate & TPE with antimicrobial & UV stabilizing additives |
| Fluid Intrusion | IP 65 - dust & water stream resistant |
| Cleaning | Stable with standard hospital surface cleaners |
| Power | 12VDC @ 3A plug-in power supply provided |
| iPad Features | Wifi, Bluetooth, cameras, speakers, microphones and screen not affected |

Regulatory

| | |
|--|---|
| Certifications (Version A) | UL/cUL 62368-1 (power supply) & FDA Registered Medical Device, Class 2 Exempt |
| Certifications (Versions B & C) | UL/cUL 1069 (nurse call) & FDA Registered Medical Device, Class 2 Exempt |
| Patents | US10863012B2 US20220066525A1 |

Accessories

| | |
|----------------------------------|---|
| Kickstand & Rail Clip | Friction hinged kickstand & springed rail clip compatible with all bed & stretcher rails as well as a wall mount |
| Security Lock | Compatible with all versions, universal key provided |

