



reddot winner

HallMonitor plus
by HATCH+MED

Epic

201A

Contact

Admit Date: 07/14/2025
 Attending: Drew Walker, MD
 Nurse: Kim H (Ph: 777-3920)

Fall Risk **Latex Allergy**

Required Before Entering Room

Visitors

- Wash Hands
- Wear Mask

Staff

- Wash Hands
- Wear Mask
- Wear Gloves
- Wear Gown

Full Code **Epic** Med Surg 2

201A

HATCH+MED

An Intelligent Display for Complex Data

Translating EMR, nurse call and bed signals into clear, actionable insights



Connecting patients to technology
in a **way that makes sense.**

Since 2012, HatchMed has been transforming hospital safety, patient engagement, and workflow efficiency by providing healthcare teams with technology designed around their real-world needs. Founded by industry expert Brian Hatch, HatchMed is built on over 30 years of combined experience in software, hardware, and design. Our agile, collaborative approach has earned the trust of prominent hospitals across the country, allowing us to adapt to diverse facility needs with impactful technology solutions.

Our mission is straightforward: to help patients leave the hospital healthier, happier, and faster. We achieve this through innovation, bringing cutting-edge solutions to patient care, while staying true to our core values of being collaborative, agile, and intuitive.

HATCH+MED

Epic Monitor curates data from the patient's record & their current stay to answer a simple question:

does the patient need your attention **right now?**

Contextual Illumination

Customizable

Patient Precautions

Configurable Chart Display

Patient Isolation Status

Assigned Care Team

Entry Requirements

201A

Contact

Admit Date: 07/14/2025
Attending: Drew Walker, MD
Nurse: Kim H (Ph: 777-3920)

Fall Risk **Latex Allergy**

Required Before Entering Room

Visitors
Wash Hands **Wear Mask**

Staff
Wash Hands **Wear Mask** **Wear Gloves** **Wear Gown**

Full Code **Epic** Med Surg 2

203A
HATCH+MED



Epic is a trademark of Epic Systems Corporation.



info@hatchmed.com

HATCH+MED



844-428-2463



Why **Epic Monitor** with HallMonitor?

Patient & Care Team Safety

Reduction in accidental intrusion
into **high isolation rooms** with
automated HallMonitor LEDs

Nurse Efficiency

Frees up hours each day for bedside
care from avoiding traditional paper
signage, **with proven ROI**

Clinical Workflow

Patient **fall risk** & RTLS **care team**
presence awareness driven by
HallMonitor



What is the cost of safety?

Our customers have helped us compare what it cost them to keep their signage as accurate as HallMonitor - **what is it costing you?**

List Price
HallMonitors

\$40 / h
Nursing Salary

\$700 / room
Cable Wiring

\$316
iPads

Real Customer Data - Contact Us for a Personalized Estimate

Savings for a 40 room med-surg unit	Time Spent per room per day	Time Spent all rooms annually	Annual Cost care team's time
Before HallMonitor Referencing EMR & manual signage updates 2x per day	11.5 minutes	2,800 hours	\$112,000
With HallMonitor Updating non-EMR statuses on screen	2 minutes	486 hours	\$19,500

HallMonitor lets your care teams focus on patients — not chores.

380 minutes / day **2,305** hours annually

Additional Patient Care Time

\$92,000 annual savings

9.8 months ROI

The Bottom Line





Staff Badge Reader

Authenticated workflows
& non-EMR field edits



Flexible Deployment

Portrait or landscape
orientation




Plug & Play Installation

Single PoE+ cable keeps the iPad charged, connected to wired internet & provides HallMonitor power. MAC whitelisting is supported.

Maintenance

Cleanable with standard medical wipes & solutions. Modular, USB C based design, built around upgradeability and seamless component replacement.



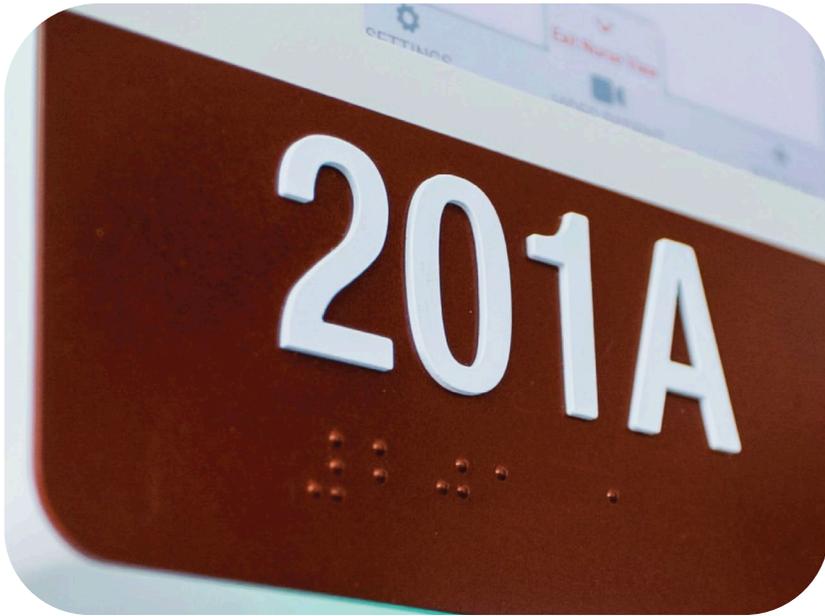

ADA Compliant & Tailor Made for You

Grade 2 Braille with raised room numbers to replace traditional signage & paper notes

Enclosure color, finish, number plate, and font — all matched to your hospital design. Software is tailored to existing workflows.

What are the customization options?

Meaningful illumination tied to real-time status — like fall risk, isolation, in-room presence, or bed safety. EMR-connected and customizable, they highlight only what matters, when it matters.

How do the LEDs impact alarm fatigue?




HallMonitor is a platform which grows with you - as iPads update, all cabling and circuitry remain constant, a simple frame upgrade is provided by HatchMed at no cost with the Upgrade Plan.

What happens when a new iPad comes out?






HallMonitor Vertical

iPad 10th & 11th (A16) Gen
HMNF-3921-KIT-PLUS

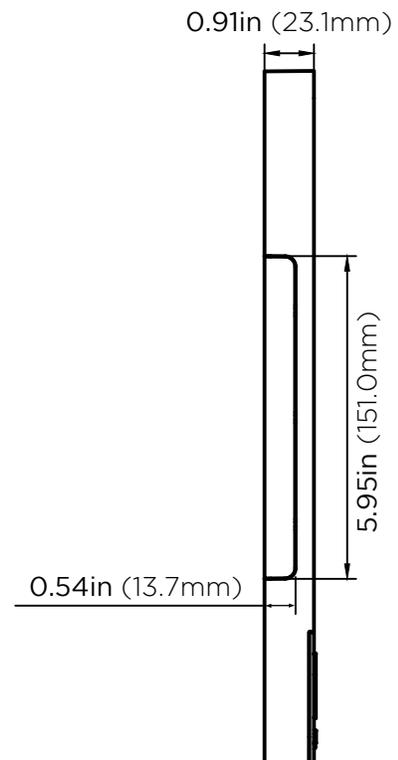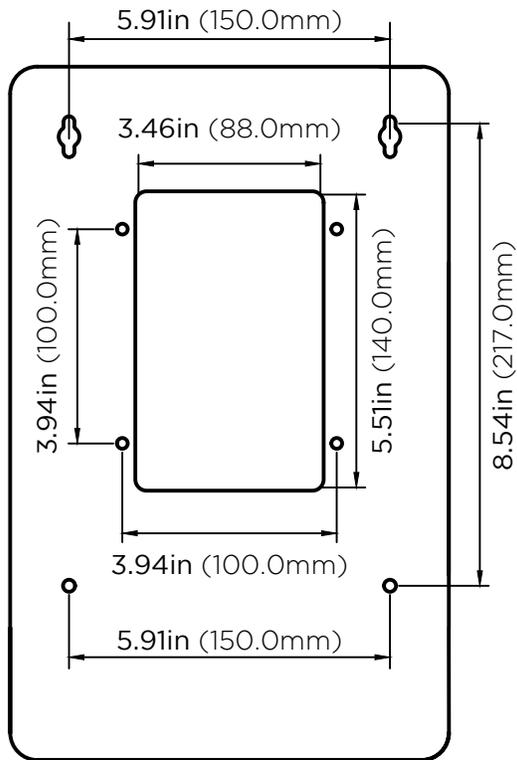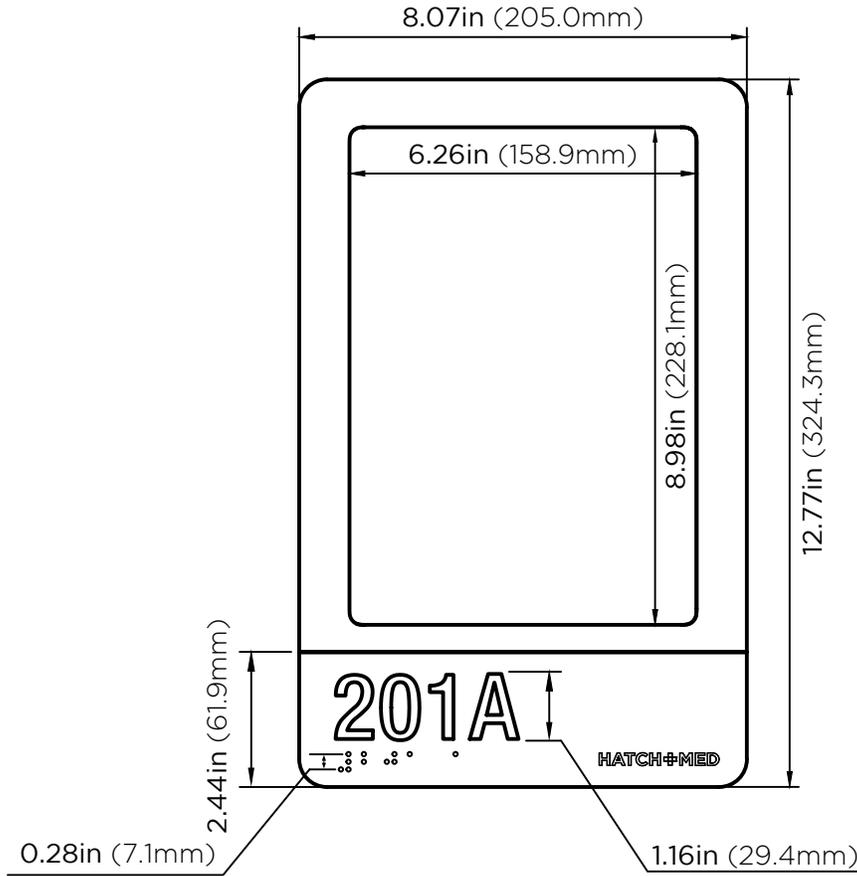

HallMonitor Horizontal

iPad 10th & 11th (A16) Gen
HMNF-3912-KIT-PLUS

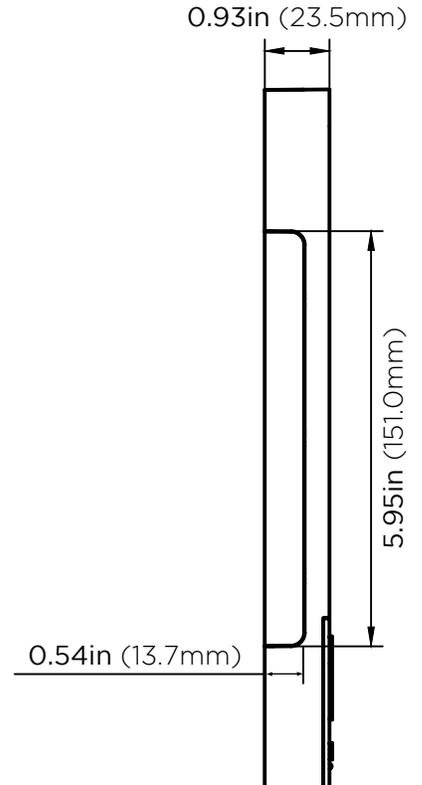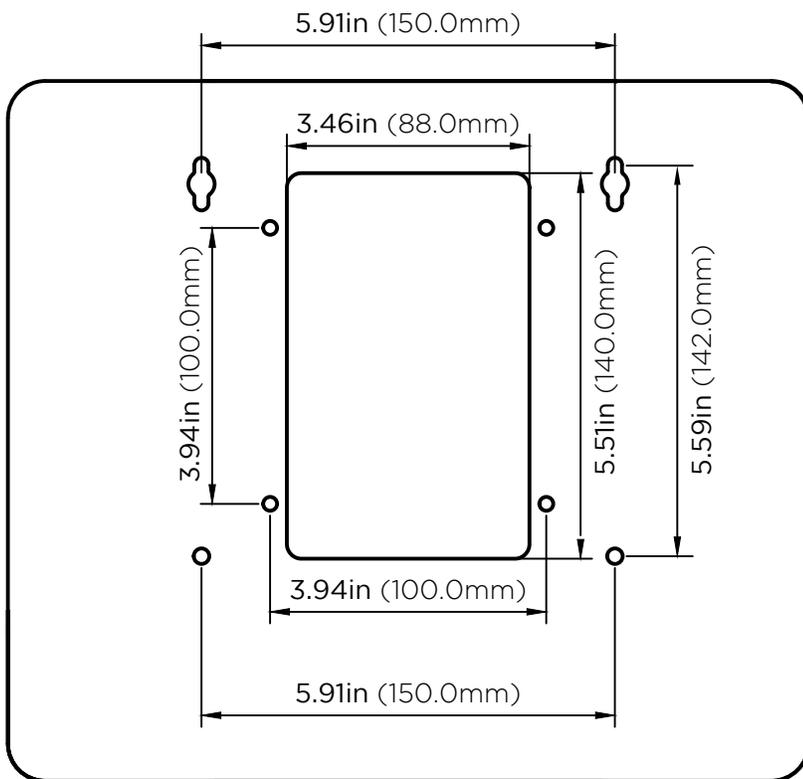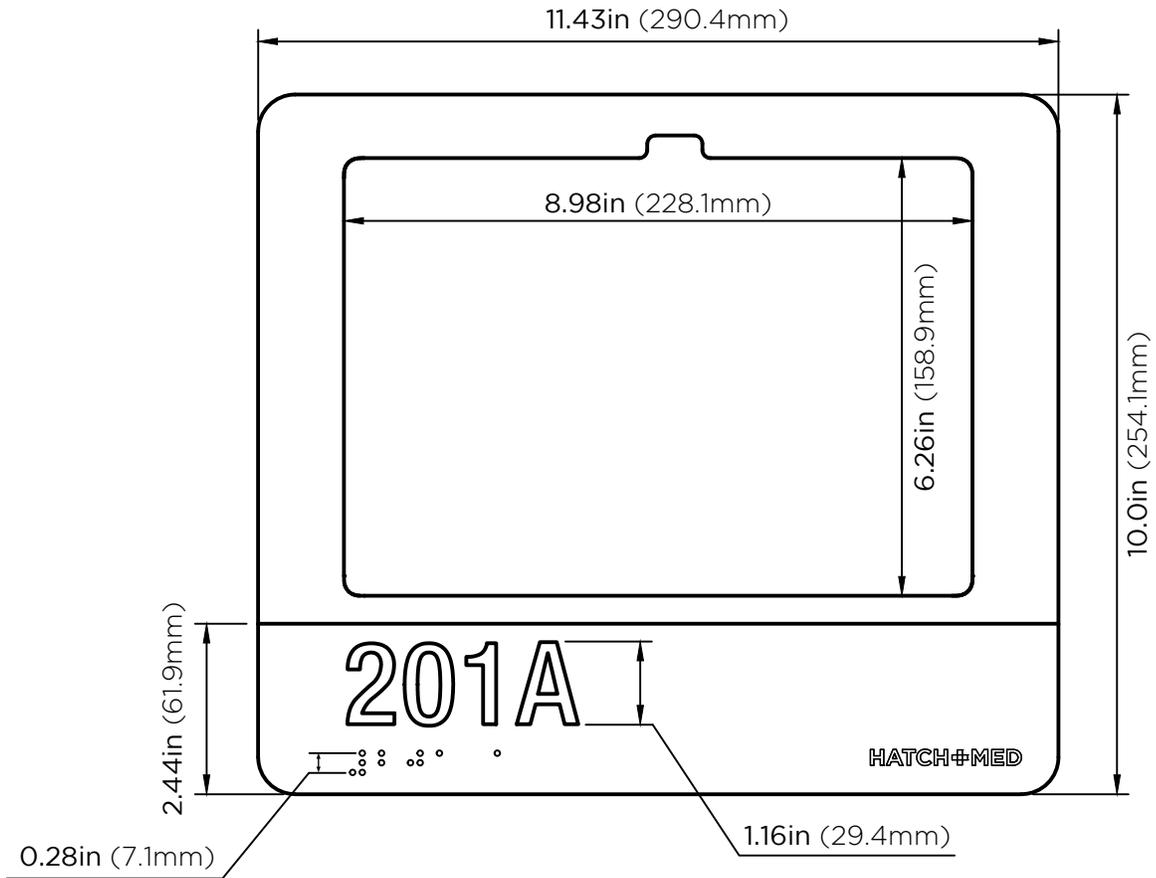

HallMonitor Specifications



Scan for info
& install guide



Epic is a trademark of Epic Systems Corporation.

Size & Weight

Devices	iPad 10th & 11th (A16) Gen
Weight	2.0 lb (without iPad) 3.1 lb (with iPad)
Dimensions	8.07" W x 12.77" H x 0.91" D (Vertical) 11.43" W x 10.00" H x 0.93" D (Horizontal)

General Info

Connection & Charging	PoE+ (Ethernet & power) or 12V @ 3A (WiFi & power)
Power Draw	Peak: ~20W (iPad fast charge) Standard operation: ~15W
Materials	UV resistant painted ABS body & Lexan illumination panels
Illumination	RGBW LED panels on all sides, support for varied brightness and patterns
RFID Badge Reader	HID OMNIKEY module supports most staff & ID badges
iPad Features	WiFi, front camera, speakers, microphones and screen not affected

Regulatory

Certifications	UL/cUL IEC 62368-1:2014 & FDA Registered Medical Device, Class 2 Exempt
Patents	US20210051223A1 US20220051532A1

Customization

Enclosure Colors	Two tone Pantone color scheme selection, custom paint matching available
Numbers & Braille	ADA compliant, 1mm raised letters and Grade 2 Braille



reddot winner



info@hatchmed.com

HATCH+MED



844-428-2463